



**FRENCH
MANUFACTURER**

maqpro



Technical Guide of Special effects

**Make-up products for simulations of wounds and various traumatisms.
Add unforgettable realism to your first aid training exercises...**



Brand's manufacturer of

Le Maquillage Professionnel – V.I.P. Make-Up – Technic Make-up – Michel Deruelle

maqpro 2 Ter rue Alasseur 75.015 Paris FRANCE

Tél : (33) 1.42.25.10.11 Fax : (33) 1.42.25.12.20 Site : www.maqpro.com e-mail : maqpro@orange.fr

S.A.S. au capital de 60.000 €- SIRET 494 440 795 00032 - RCS PARIS 494 440 795

APE 2042 Z - N° TVA Intracommunautaire FR 46 494 440 795



First aid trauma & casualty simulations kits

Le Maquillage
Professionnel

2 models suggested

Ref. KM1

Including 23 products

- | | |
|--|-----------------------------------|
| 1 DVD-ROM – Teaching aid support | 1 COLLODION OFFICINAL 10ml |
| 1 PLASTO-NAT 29g | 1 LATEX 30ml |
| 1 NORMAL MAKE-UP BLOOD 125ml | 1 SEALOR 10ml |
| 2 BRUISES SPONGES | 1 MAKE-UP BLOOD JELLY Tube 30g |
| 2 FOAM LATEX SPONGES | 1 ADHESIVE AND LATEX SOLVENT 30ml |
| 1 BLISTERS SIMULATION PRODUCT Tube 30g | 1 JELLY MAKE-UP REMOVER 30ml |
| 3 CREAM BLUSHERS in trio limps yellow/blue/blood red | 1 WOOD SPATULA |
| 3 CREAM BLUSHERS in trio limps paleness/red/khaki | 1 FOAM SPONGE <i>round shape</i> |
| 1 PERSPIRATION LIQUID 30 ml | |



Rigid Case with foam of protection

Ref. KM2

Including 28 products

- | | |
|--|-----------------------------------|
| 1 DVDROM – Teaching aid support | 1 SEALOR 10ml |
| 1 PLASTO-NAT 29g | 1 MAKE-UP BLOOD JELLY Tube 30g |
| 1 NORMAL MAKE-UP BLOOD 125ml | 1 ADHESIVE AND LATEX SOLVENT 30ml |
| 2 BRUISES SPONGES | 1 MAKE-UP REMOVER JELLY 30ml |
| 2 FOAM LATEX SPONGES | 1 WOOD SPATULA |
| 1 BLISTERS SIMULATION PRODUCT Tube 30g | 1 FOAM SPONGE <i>round shape</i> |
| 3 CREAM BLUSHERS in trio limps yellow/blue/blood red | 1 SPATULA stainless |
| 3 CREAM BLUSHERS in trio limps paleness/red/khaki | 1 BRUSH MAKE UP N°2 |
| 1 PERSPIRATION LIQUID 30 ml | 1 ERCEFILM |
| 1 COLLODION OFFICINAL 10ml | 1 PLASTO NAT 17ml |
| 1 LATEX 30ml | 1 LATEX PROTECTIVE COVERING 10 ml |



Rigid case with foam of protection



FREE DVD-ROM – TEACHING AID SUPPORT COMPRISING:

- technical resumes of 25 simulation cases
- videos
- more than 100 mo of technical descriptions, advices & tips on how to use the products





Réf. **KM3**

Including 21 products

- | | |
|----------------------------------|----------------------------------|
| 1 DVD-ROM – Teaching aid support | 1 BLUE CREAM BLUSHER 7ml |
| 3 PLASTO-NAT (29gx3) | 1 BLOOD RED CREAM BLUSHER 7ml |
| 1 THICK MAKE-UP BLOOD 250ml | 1 FOAM SPONGE <i>round shape</i> |
| 1 FLUID MAKE-UP BLOOD 250ml | 1 BRUISES SPONGE |
| 1 MAKE-UP BLOOD JELLY Tube 30 g | 1 WOOD SPATULA |
| 1 PALENESS CREAM BLUSHER 15ml | 2 WOOD STICKS |
| 1 PERSPIRATION LIQUID 60ml | 1 PAINTBRUSH (pony hair) |
| 1 LATEX 60ml | 1 SEALOR 15ml |
| 1 YELLOW CREAM BLUSHER 7ml | 1 LATEX PROTECTIVE COVERING 15ml |



*Rigid case with foam
of protection*

This kit is perfect for exercises and operations.

Réf. **KM4**

DESCRIPTION of the kit

Including 32 products

**For this kit we propose you 2
different packaging:**



*a fabric vanity-case
Réf. KM4 / A*

*or Rigid case with foam
of protection
Réf. KM4 / B*

- 1 DVD-ROM – Teaching aid support
- 5 PLASTO-NAT (29gx5)
- 1 THICK MAKE-UP BLOOD 250ml
- 1 FLUID MAKE-UP BLOOD 250ml
- 1 COAGULATED BLOOD Tube 30g
- 1 MAKE UP BLOOD JELLY Tube 30g
- 1 PALENESS CREAM BLUSHER 15ml
- 1 PERSPIRATION LIQUID 60ml
- 2 BLISTERS SIMULATION PRODUCT Tube 30g
- 1 LATEX 60ml
- 1 YELLOW CREAM BLUSHER 7ml
- 1 BLACK CREAM BLUSHER 7ml
- 1 BLOOD RED CREAM BLUSHER 7ml
- 1 BLUE CREAM BLUSHER 7ml
- 1 RED CREAM BLUSHER 7ml
- 2 FOAM SPONGES *round shape*
- 2 BRUISES SPONGES
- 1 WOOD SPATULA
- 2 WOOD STICKS
- 1 PAINTBRUSH (pony hair)
- 1 MAKE-UP REMOVER MILK 250ml
- 1 SOILING POWDER 50ml
- 1 SEALOR 15ml
- 1 LATEX PROTECTIVE COVERING 15ml



**FREE DVD-ROM – TEACHING AID SUPPORT
COMPRISING:**

- **technical resumes of 25 simulations cases**
- **videos**
- **more than 100 mo of technical descriptions,
advices & tips on how to use the products**





MSP 01
Ref. MSP 01

Items	Cont.	Quant.
DVD-ROM – Teaching aid support		1
PALETTE of 10 CREAM BLUSHERS *	50 ml	1
PALNESS CREAM BLUSHER	50 g	1
NORMAL MAKE-UP BLOOD	60 ml	1
FLUID MAKE-UP BLOOD	60 ml	1
THICK MAKE-UP BLOOD	60 ml	1
MAKE-UP BLOOD JELLY (ointment)	30 g	1
LATEX PROTECTIVE COVERING	60 ml	1
LATEX	30 ml	1
COLLODION OFFICIAL	30 ml	1
PERSPIRATION LIQUID	30 ml	1
MASTIC SPIRIT GUM	30 ml	1
PLASTO-NAT	85 ml	1
ADHESIVE AND LATEX SOLVENT	30 ml	1
MAKE-UP REMOVER JELLY	200 ml	1
FOAM SPONGE		2
BRUISES FOAM SPONGE	Bath of 3	1
SEALOR	30 ml	1

*Cream blushers colours: black-white-yellow-red-blue-blood red-khaki-pale flesh tone-medium flesh tone-dark flesh tone



MSP 02
Ref. MSP 02

Item	Cont.	Quant.
DVD-ROM – Teaching aid support		1
PALETTE of 10 CREAM BLUSHERS *	50 ml	1
PALNESS CREAM BLUSHER	50 g	1
NORMAL MAKE-UP BLOOD	60 ml	1
FLUID MAKE-UP BLOOD	60 ml	1
THICK MAKE-UP BLOOD	60 ml	1
MAKE-UP BLOOD JELLY Tube	30 g	1
LATEX PROTECTIVE COVERING	60 ml	1
LATEX	30 ml	1
COLLODION OFFICIAL	30 ml	1
PERSPIRATION LIQUID	30 ml	1
MASTIC SPIRIT GUM	30 ml	1
PLASTO-NAT	85 ml	1
ADHESIVE AND LATEX SOLVENT	30 ml	1
MAKE-UP REMOVER JELLY	200 ml	1
FOAM SPONGE		2
BRUISES FOAM SPONGE	Bath of 3	1
SEALOR	30 ml	1
SPATULA STAINLESS		1
PAINTBRUSH N° 2		1
FIXATING POWDER	50 ml	1
BLISTER SIMULATION PRODUCT Tube	30 g	1
BIOLOGICAL CLEANER	60 ml	1

*Cream blushers colours: black-white-yellow-red-blue-blood red-khaki-pale flesh tone-medium flesh tone-dark flesh tone



MSP 03

Ref. MSP 03

Désignation	Cont.	Quant.
DVD-ROM – Teaching aid support		1
PALETTE of 10 CREAM BLUSHERS *	50 ml	1
PALNESS CREAM BLUSHER	50 g	1
NORMAL MAKE-UP BLOOD	60 ml	1
FLUID MAKE-UP BLOOD	60 ml	1
THICK MAKE-UP BLOOD	60 ml	1
MAKE-UP BLOOD JELLY Tube	30 g	1
LATEX PROTECTIVE COVERING	60 ml	1
LATEX	30 ml	1
COLLODION OFFICINAL	30 ml	1
PERSPIRATION LIQUID	30 ml	1
MASTIC SPIRIT GUM	30 ml	1
PLASTO-NAT	85 ml	1
ADHESIVE AND LATEX SOLVENT	30 ml	1
MAKE-UP REMOVER JELLY	200 ml	1
FOAM SPONGE		2
BRUISES FOAM SPONGE	Bath of 3	1
SEALOR	30 ml	1
SPATULA STAINLESS		1
MAKE UP BRUSH N° 2		1
FIXATING POWDER	50 ml	1
BLISTER SIMULATION PRODUCT Tube	30 g	1
CLEANER	60 ml	1
PAINTBRUSH N° 6		1
SIMULATION SKIN PRODUCT Tube	30 g	1
COAGULATED BLOOD Tube	30 g	1

*Cream blushers colours: black-white-yellow-red-blue-blood red-khaki-pale flesh tone-medium flesh tone-dark flesh tone



KIT ECOLE N° 1

Ref. KIT 01

Désignation	Cont.	Quant.
DVD-ROM – Teaching aid support		1
PALETTE of 10 CREAM BLUSHERS *	50 ml	1
PALNESS CREAM BLUSHER	50 g	1
NORMAL MAKE-UP BLOOD	1 000 ml	1
FLUID MAKE-UP BLOOD	1 000 ml	1
THICK MAKE-UP BLOOD	1 000 ml	1
MAKE-UP BLOOD JELLY Tube	90 g	1
LATEX PROTECTIVE COVERING	125 ml	1
LATEX	125 ml	1
COLLODION OFFICINAL	125 ml	1
PERSPIRATION LIQUID	125 ml	1
MASTIC SPIRIT GUM	60 ml	1
PLASTO-NAT	200 ml	1
ADHESIVE AND LATEX SOLVENT	500 ml	1
MAKE-UP REMOVER JELLY	200 ml	1
FOAM SPONGE		4
BRUISES FOAM SPONGE	Bath of 3	2
SEALOR	125 ml	1
SPATULA STAINLESS		1
MAKE UP BRUSH N° 2		1
FIXATING POWDER	45 g	1
BLISTER SIMULATION PRODUCT Tube	90 g	1
CLEANER	500 ml	1
PAINTBRUSH N° 6		1
SIMULATION SKIN PRODUCT Tube	90 g	1
ERCEFILM	150 ml	1

*Cream blushers colours: black-white-yellow-red-blue-blood red-khaki-pale flesh tone-medium flesh tone-dark flesh tone



FAIN FEELING

Made up with

- Cleaner (p.12)
- Paleness cream blusher (p.9)
- Sponge (page 9)
- Khaki cream blusher (p.9)
- Bruises sponge (p.9)
- Perspiration liquid (p.9)

Make-up remover

- Make-up remover jelly or Make-up remover milk (p.12)
- Cleaner (p.12)

CYANOSIS

Made up with

- Cleaner (p.12)
- Paleness cream blusher (p.9)
- Sponge (p.9)
- Khaki cream blusher (p.9)
- blue cream blusher (p.9)

Make-up remover

- Make-up remover jelly or Make-up remover milk (p.12)
- Cleaner (p.12)



HEAVY BURN

Made up with

- Latex protective covering (page 8)
- Red cream blusher (page 9)
- Blood red cream blusher (page 9)
- Yellow cream blusher (page 9)
- Khaki cream blusher (page 9)
- Black cream blusher (page 9)
- Blister simulation product (page 8)
- ERCEFILM (page 10)
- Artificial blood (page 11)

Make-up remover

- Adhesive and latex solvent (p.12)
- Cleaner (p.12)

THERMAL BURN

Made up with

- Latex protective covering (p.8)
- Red cream blusher (p.9)
- Latex (p.10)
- Black cream blusher (p.9)
- Blister simulation product (p.8)
- Make-up blood (p.11)

Make-up remover

- Adhesive and latex solvent (p.12)
- Cleaner (p.12)

BASIC BURN

Made up with

- Red cream blusher (p.9)
- Blood red cream blusher (p.9)
- Yellow cream blusher (p.9)
- Collodion Officinal (p.10)

Make-up remover

- Adhesive and latex solvent (p.12)
- Cleaner (p.12)





basic wounds

Simulation of commonly found wounds-flesh wound-open wound- glass embedded wound-cut finger...

Made with

- Latex protective covering (p.8)
- Plasto-Nat (p.8)
- Sealor (p.8)
- or
- ERCEFILM (p.10)
- Cream blushers :
yellow – blue – blood red (p.9)
- Make-up blood jelly (p.11)
- Make-up blood (p.11)

Make-up remover

- Adhesive and latex Solvent (p.12)
- Cleaner (page 12)

Serious wounds

Scalp, eviscerations, throat cut...

Made with

- Latex protective covering (p.8)
- Skin simulation product(p.10)
- Latex (p 10)
- Powder (tarif Réf. 1062)
- Blood red cream blusher (p 9)
- Make-up blood (p 11)

Make-up remover

- Adhesive and latex Solvent (p.12)
- Cleaner (p.12)



MANŒUVRES ET PLAN ROUGE

To create large wounds on multiple victims of disaster for instance we recommend using prosthesis:



Ref : PR01



Ref : PR05



Ref : PR04



Ref : PR02



Ref : PR06



Ref : PR07



Ref : PR03



Ref : PR08



PLASTO-NAT

The PLASTO-NAT is a modelling paste used like modelling clay on the skin (use the LATEX PROTECTIVE COVERING ref.2070 to 2073 before the application of it). Depending of the effect wanted you can make it stronger adding a thin couch of SEALOR which is allowed a compression without to put out of shame your realization.

Once the volume puts in place you can make up the PLASTO-NAT with the CREAM BLUSHERS.

➔ *Make-up remover:* We recommend to use the MAKE-UP REMOVER JELLY (ref.1041 to 1043/A) to get rid of the colouring and to take off everything use a cotton wool and adhesive and latex Solvent (réf.2023 to 2023 /C).



Ref. 2026 17ml



Ref. 2027 29 Grsl



Ref. 2027/A 85ml



Ref. 2027/B 200ml

➔ *Les couleurs :*



Pale flesh



Medium flesh



Brown flesh

SEALOR

This is a transparency film which allows making the PLASTO-NAT stronger. We recommend spreading it with a paintbrush or the fingers.

Once the wound becoming strong it can resist to manual compressions, water and you can make it up with make-up blood or spread some perspiration liquid...

It doesn't make any skin irritation.

➔ *Make-up remover:* With a cotton wool and adhesive and latex Solvent (réf.2023 to 2023 /C), use it also to clean the paintbrush.

Caution: Never apply near the eyes or the mucous membranes.



Réf. 2035 60ml



Réf. 2036 125ml



Réf. 2035/B 500ml



Réf. 2036/B 1 000ml

LATEX PROTECTIVE COVERING

We recommend applying this protective on the skin before to create any realization, except "feeling paint"- cyanosis - plain burn and close to the eyes, in order to protect the skin from the allergic reactions. Leaves the skin breathing and easily removes without pull out the hair.

➔ *Easy way:* It makes the PLASTO-NAT lasting longer

➔ *Make-up remover:* With a cotton wool and adhesive and latex Solvent (réf.2023 to 2023 /C).



Réf. 2070 60ml



Réf. 2071 125ml



Réf.2072 500ml



Réf. 2073 1 000ml

Caution: Never apply near the eyes or the mucous membranes.

BLISTER SIMULATION PRODUCT

This item allows simulating burns with closed blisters, after having apply the LATEX PROTECTIVE COVERING (Ref.2070 to 2073) you can colour the BLISTER SIMULATION PRODUCT with the CREAM BLUSHERS (ref.1000 to 1004), create the blisters and fix them with the RCFILM (ref.RC001) or with the LATEX (ref.2009 to 2012).

➔ *Make-up remover:* With a cotton wool and adhesive and latex Solvent (réf.2023 to 2023 /C).



Réf. 2150 30 G



Réf. 2151 90 G



Réf. 2052 500ml



CREAM BLUSHER

The CREAM BLUSHER is the standard product used by the professionals of the make-up for all their achievements.

The CREAM BLUSHER resists water and "holds" 12 hours without problem - does not stain clothing.

It will be used to make up the achievements of your wounds out of Latex or Plasto-Nat, for the effects of volumes and the connections with the skin.

COLORS



The cups can be gathered by 10 in a plastic pallet (réf. 1003)

* Composition of the pallet cream blusher: *black-white-yellow-red-blue-blood red-khaki-pale flesh tone-medium flesh tone-dark flesh tone*
You can also have this product in plate of 50g (réf. 1000) for the make-up of great surfaces (face, body) ideal for make-up paleness.

➔ *Make up remover*: with the **make up remover jelly** Réf. 1041, the **Cleaner** Réf. 1034 or the **make up remover milk** Réf. 1038 /A



Réf. 1000
Li of 50g



Réf. 1003 50ml pallet OF 10 colors *
Cups Réf. 1003.A 5ml
reload for pallet



Réf. 1004
17ml



Réf. 1004 / A
7ml



Réf. 1003 / C 6 ml
N° 1 paleness/red/khaki N° 2 yellow/blue/blood-red N° 3 yellow/blue/black N° 4 red/paleness/black



SPONGES

Essential to apply the CREAM BLUSHER when you are not use a brush. The triangular form allows going in the hollows while the round form allows the make-up of a wider surface. The choice of the shape of sponge also depends on the "turn of hand" of each one.

After each use, wash the sponge with water and a little household soap.



Réf. 1054



Réf. 1055/A Réf. 1055/B

BUISES SPONGE

With the bruisse sponge it makes possible to realise the effect of vessels burst on the point of impact, for this effect, one uses the CREAM BLUSHER can also be used for to make beads of sweat with the perspiration liquid.



Réf. 1061 / A

PERSPIRATION LIQUID

Liquid used to simulate perspiration. Use it with the bruisse sponge (réf. 1061) to imitate the beads of sweat.

➔ *Make-up remover*:
Remove it with water.



500ml 125ml 60ml
Réf. 2029/D 2029 2029/A



OTHERS

ERCEFILM - LATEX - COLLODION OFFICINAL - SIMULATION SKIN PRODUCT - MASTIC SPIRIT GUM

LATEX

It acts of natural Latex pre-vulcanized which is worked very easily, while drying it becomes brilliant and/or transparent. It uses to make simulation of burns or the effects of scraps of skin etc.... You can colour latex with CREAM BLUSHER.

➔ *Easy way:* Apply the latex protective covering before the latex to limit the risks of allergic reactions (see page N°8).

➔ *Make-up remover:* With a cotton wool and adhesive and latex solvent (réf.2023 to 2023 /C).

Caution: Never apply near the eyes or the mucous membranes.



1 000ml Réf. 2012
500ml Réf. 2011
125ml Réf. 2010
60ml Réf. 2009

MASTIC SPIRIT GUM

This adhesive makes it possible to fix all the false hairs.

➔ *Make-up remover:* ➔ *Make-up remover:* With a cotton wool and adhesive and latex solvent (réf.2023 to 2023 /C) in order to not leave of red trace on the skin and to not drain it.

Caution: Never apply near the eyes or the mucous membranes.



60ml Réf. 2016
30ml Réf. 2015

COLLODION OFFICINAL

It uses to simulate scars or burns. It will draw the skin while drying to form according to the application a more or less major effect. As it is transparent, you can add colour on the skin with the make-up cream to give an effect of depth.

➔ *Make-up remover:* With a cotton wool and adhesive and latex solvent (réf.2023 to 2023 /C).



60ml Réf. 2024
125ml Réf. 2024/E
500ml Réf.2024/D
1 000ml Réf. 2024/B

Caution: Never apply near the eyes or the mucous membranes.

SIMULATION SKIN PRODUCT

This product can simulate « human skin ». Useful to realize: an abdomen wound, a sliced throat, etc... This effect of “skin” with latex can be sewn. For the coloration, you can use the make-up creams.

➔ *Make-up remover:* With a cotton wool and adhesive and latex solvent (réf.2023 to 2023 /C).



Réf. 2180 30 Gr

Réf. 2181 90 Gr

Réf. 2082 500ml


ERCEFILM

Pellicular bandage in spray.

This product replaces latex for the fixing of the blisters simulation product.

➔ *Make-up remover:* With a cotton wool and adhesive and latex solvent (réf.2023 to 2023 /C).

Caution: Never apply near the eyes or the mucous membranes.

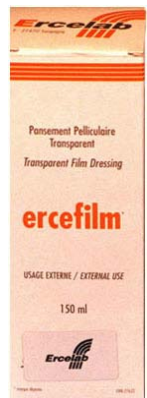

150ml Réf. RC001



MAKE-UP BLOOD (bottle)

The **MAKE-UP BLOOD** it's an exclusive creation of our manufacture and brings lot of valuable qualities:

- a very realistic colour of the real blood.
- it's a food product, it means you can make it in contact with the skin and even put it in the mouth without any danger
- it doesn't stained
- it can be easily cleaned with some cold water and soap on the skin and the cloth.

The MAKE-UP BLOOD BOTTLE exists in more or less thick versions:
L= Liquid N= Normal E= Thickly

**DON'T DILUTE THE MAKE-UP BLOOD
OR IT WILL LOSE ITS ATTRIBUTES**




5 000ml Ref. 2005/A

COAGULATE BLOOD

Pasty version of artificial blood intended to simulate scratches or coagulated wounds.
It dries quickly and holds well.

➔ *Make up remover:* Cleaned with tepid water or the CLEANER (Ref.1034 to 1037/A).

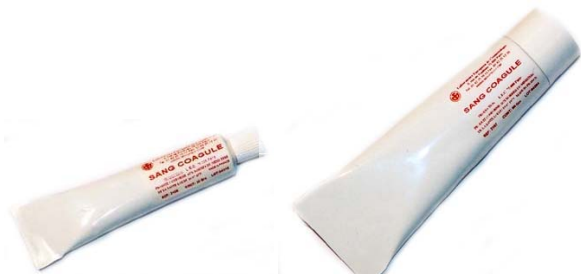

Ref. 2106 30 G

Ref. 2107 90 G



Ref. 2007/A 500ml

MAKE UP BLOOD JELLY

This texture makes a great effect of flesh with sharp; it doesn't make filament like the COAGULATED BLOOD and doesn't stick; it can be worked in transparency or thickness and create for instance a cut finger, a wound with flesh with sharp, a cut throat...

➔ *Make up remover:* Cleaned with tepid water or the CLEANER (Ref.1034 to 1037/A).



Ref. 2130 30 G

Ref. 2131 90 G



Ref. 2031/A 500ml



MAKE-UP REMOVER

Le **Maquillage**
Professionnel

Cleaner

It's a very soft make-up remover. alcohol, detergent and perfume free, which allows you to get ride of all the coloration made with the CREAM BLUSHERS.

It lets the skin clear and not greasy.

Also allows to withdraw the make-up tasks on clothing

Réf.	Contenance
1034	Plastic bottle 60 ml
1035	Plastic bottle 125 ml
1036	Plastic bottle 500 ml
1037	Plastic bottle 1 000 ml
1037/A	Plastic bottle 5 000 ml

Make-up remover milk

Our MAKE-UP REMOVER MILK is a very soft product which can be used on all kind of skins. It leaves a light greasy cover to moisturize the skin.

Réf.	Contenance
1038/A	Plastic bottle 60 ml
1038	Plastic bottle 125 ml
1039	Plastic bottle 500 ml
1040	Plastic bottle 1 000 ml

Make-up remover jelly

- Recommended for sensible skins;
- To use directly on the skin and then wipe gently with a tissue;
- Leaves a light greasy cover that you can keep it to moisturize the skin or clean it with the CLEANER or some water and soap.

Réf.	Contenance
1041	Plastic jar 200 ml
1043	Plastic jar 500 ml
1043 / A	Plastic jar 5 000 ml
1043 / B	Plastic jar 10 000 ml

Ahesive and latex solvent

This is the essential product to have as it allows you to clean up with wool cotton any trace of PLASTO-NAT, LATEX, COLLODION, SEALOR... It leaves soft skin without redness or dryness. The ADHESIVE AND LATEX SOLVENT can remove all the glues.

Caution: Never apply near the eyes or the mucous membranes.



Paintbrushes cleaner

It's a make-up remover specially made to clean gently the paintbrushes and keep them longer. Just put the paintbrush in the liquid for few seconds then dry it softly with a tissue and then your paintbrush is clean, disinfected, dry and ready to use.

Réf.	Contenance
2037	Glass bottle 60 ml
2038	Glass bottle 125 ml
2039	Can 500 ml
2040	Can 1 000 ml